

putting ON PAPR

Donning Checklist 1/3

PAPR system

coveralls

two pair long
cuff gloves


battery
pack

belt

gown

shoe
liners

booties

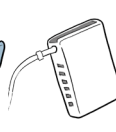

White-hooded PAPR Series 2080

Shoe liners + gloves

1 Qualified monitor is present.

2 Set up cart for donning and doffing.

3 Remove personal effects (phone, keys, badges, pens, etc.).

4 Put on clear inner shoe liners over your shoes and pants.

5 Put on extended cuff inner gloves.

Battery + Helmet

6 Secure belt snug around waist with clip opening upwards. Place excess belt into clip and loop.

7 Secure the battery on the left side of the belt.

8 Connect helmet to battery.

9 Check that 3 LED lights are green to ensure battery is fully charged.

10 Adjust temple holder (top hole for larger head, 2nd hole works for most though).

11 Put helmet on head and size to fit by turning ratchet knob to adjust the headband. Ensure long hair is in a low bun or ponytail is above shoulders.

12 Bend over to ensure helmet is tight enough to stay on without slipping.

13 Adjust fan speed to comfort.

14 Disconnect the battery and remove helmet.

PAPR Donning Checklist 2/3

Hood

15

Open the hood bag.



16

Set aside outer cap.



17

Orient the hood so markings face you and are in upward direction and flat edge is at top.



18

Hold the upper corners and **shake** the hood until it fully unfolds.




19

Connect the center hole of the hood to the center hook on the front of the helmet. Once you are hooked, **keep pushing** until you hear a snap.



20

Connect the right and left sides of the hood to the helmet in the same way until you hear a snap.

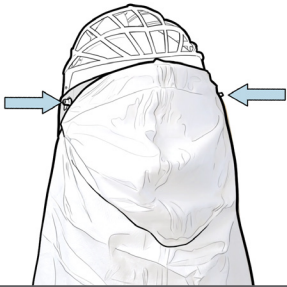

21

Grab the bottom of the hood and **pull** it back and over the helmet until Velcro is visible.



22

Take the cap and **fit it** fully over the top of the hood.



23

Set assembled PAPR hood aside.



Coveralls, Tape, Booties

24

Put on coveralls outside room. Feet in first. Open-end faces your front. **Do not** zip up yet.



25

Ensure proper mobility of arms and legs within coveralls.



26

Tape coverall and wrist 50/50 with 1.5 - 2 ft of masking tape. And form a release tab at the end of the loop.



27

Don booties.




PAPR Donning Checklist 3/3

Zip and Tie

28 **Sanitize** gloves and touched dispenser.




29 **Insert** PAPR power cord connector into battery.

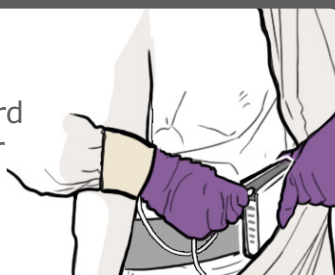

30 **Put** the helmet on, letting the hood drape around your shoulders over the coveralls.



31 **Insert** fingers into facial liner and **tug** it around your face until it tucks under your chin. Last chance to **adjust** the helmet!

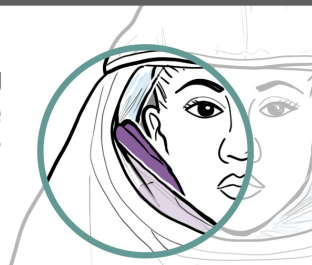

32 **Tuck** cord into your suit.




33 **Zip up** the front of the coveralls.




34 **Find** both neck ties and bring them to the front and **tie** a shoelace knot close to neck.



35 **Find** both torso ties from center back of hood, bring them around and **insert** them through the front holes.

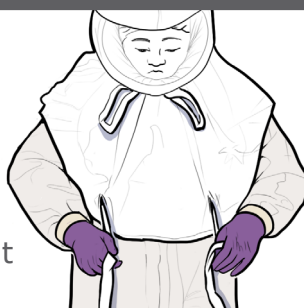

36 **Tie** into shoelace knot.

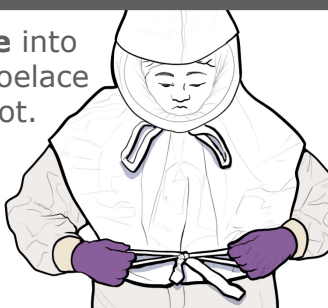

Gown

37 **Unfold** and **put on** gown.



38 **Secure** neck of gown, making sure neck tie of hood is covered.




39 **Tie** gown at waist at left side.



40 **Put on** extended cuff outer gloves, **tuck** gown sleeves into gloves.



41 Safety monitor completes **360° check** of PPE.

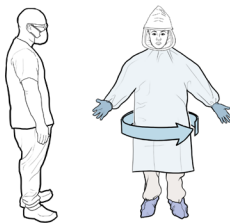


Entry

42 **Sanitize** gloves and touched dispenser.




43 **ENTER room** **DO NOT** enter the room without a proper fit of the PAPR system.

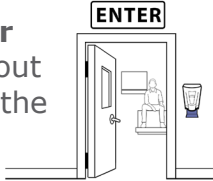





44 **Do not** touch or re-adjust **PAPR** system inside the room.



Setup, wipe

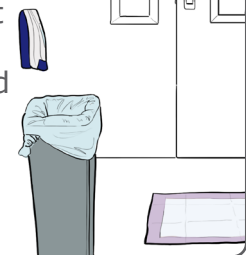
- 1 **Make sure** that you have everything you need to complete a successful doffing. (Trash can, hand sanitizer, balance aid). 
- 2 **Check** that your partner outside of the room is aware you are preparing to doff. 
- 3 **Inspect** PPE for visible signs of contamination. 
- 4 Have your **partner inspect** your back through the window on the door. 
- 5 **Wipe away** contamination with SaniWipes. 
- 6 **Place** a large underpad in front of the door. 
- 7 **Ensure** that you have sanitizer and trash can within reach of pad. 
- 8 **Sanitize** gloves and touched dispenser. 

Cap + Gown

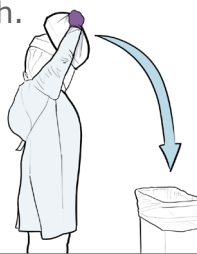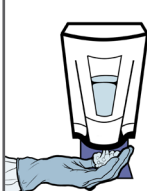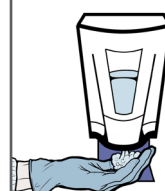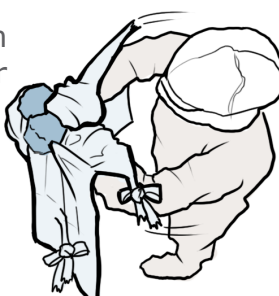- 9 With one hand, **reach** up to the back of the helmet, **pinch** and **remove** cap from top of PAPR hood, throw into trash. 
- 10 **Sanitize** gloves and if touched, dispenser. 
- 11 Gently **pull** the sleeves of the gown out from under each outer glove cuff, working your way around the wrist. 
- 12 **Sanitize** gloves and if touched, dispenser. 
- 13 **Grasp** gown at sides over waist tie, **pull forward** to break tie. 
- 14 **Cross** your arms in an "X" to **grasp** gown at upper arm, close to shoulder, **avoid** touching PAPR hood. 
- 15 **Uncross** your arms slowly, extending them away from your body, inverting the gown. **Extend** your arms outward until they are straight. 
- 16 Keeping arms extended away from your body, **bundle** the gown and discard. 

Booties


17 **Sanitize** gloves and if touched, dispenser.

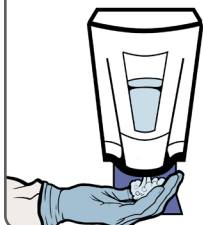

18 Using balance aid if appropriate, with one hand **remove** one bootie. As you remove it, step with that foot onto pad. Throw bootie away.




19 **Remove** other bootie. Step onto pad. Throw bootie away.

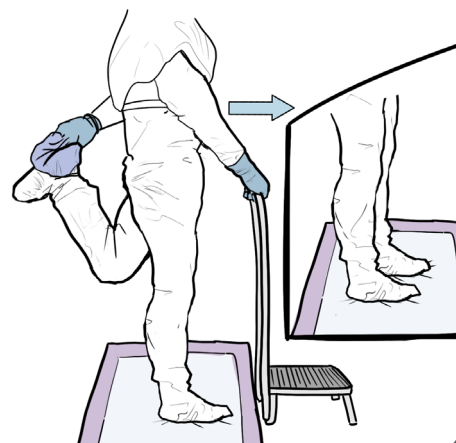


Gloves


20 **Sanitize** gloves and if touched, dispenser.

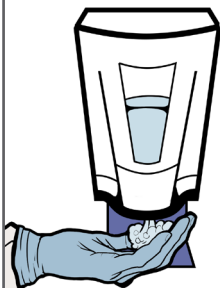

21 **Pinch** glove cuff and pull over fingers, but do not remove glove. With covered fingers, **pinch** and **pull** opposite glove off. Using other hand, **remove** initial glove and discard both gloves.

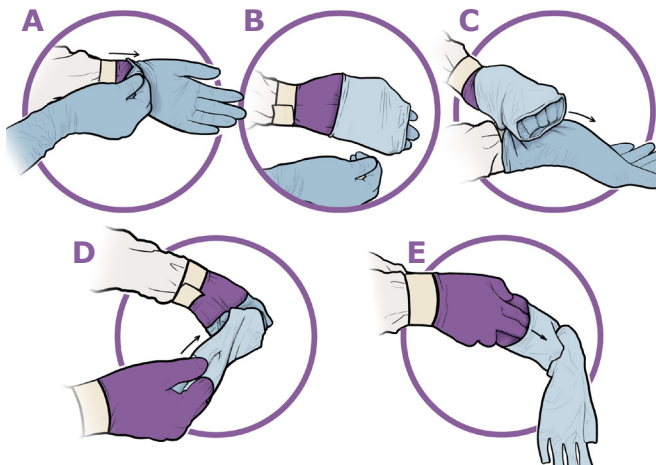

22 **Sanitize** gloves and if touched, dispenser.

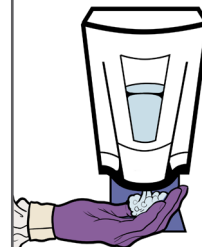


Ties

23 **Untie** torso ties, **ball up** and **pull through** one side, and discard.




24 **Sanitize** gloves and if touched, dispenser.

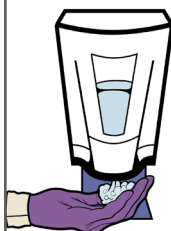

25 **Untie** neck ties, **ball up** and **pull through** one side, and discard.



26 **Sanitize** gloves and if touched, dispenser.

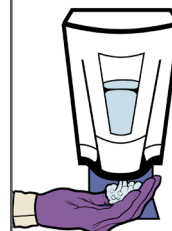

Coveralls


27

In a controlled manner, **pull** tape down and away to remove tape from both wrists.



28

Sanitize gloves and if touched, dispenser.

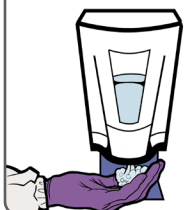

29

Avoiding touching the inside of the PAPR hood, **zip down** coveralls with one hand.




30

Grasp the coveralls near lower back, at the hips, and **pull** the coveralls straight down.



31

Touching only the outside of the coveralls, **push** the coveralls down below the knees.

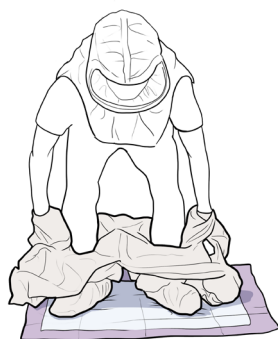

32

Step on the crotch of the suit to remove feet from coverall. **Repeat** step motion until feet are free.




33

Step off coveralls and leave the coveralls for the next person to pick up.



34

Sanitize gloves and if touched, dispenser.

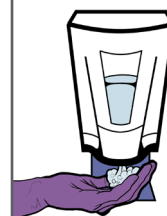

Exit

35

Observer opens door.



36

Exit patient room.



Hood

37 **Sanitize** gloves and touched dispenser.




38 **Place** both hands under the front of the hood, palms forward. **Do not cross arms.** Reach to the back of the helmet on each side. Fully **gather** hood material and **tightly roll** touching only the inside of the hood.



39 Continue to **roll** the hood material up across the back and top of the helmet until you reach the top of the head, in line with the ears.



40 **Move** fingers to the edge of the face shield and unsnap the temple post from the helmet and fold visor in half.

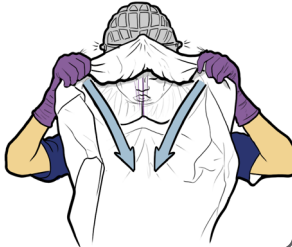


41 With an upward motion, **pull** folded visor from center post.



Remove hood.


DO NOT TOUCH YOUR FACE!

42 **Sanitize** gloves and if touched, dispenser.

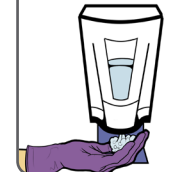


Gloves + Observer

43 **Pinch** glove cuff and **pull** over fingers, but **do not remove** glove. With covered fingers, **pinch** and **pull** opposite glove off. Using other hand, **remove** initial glove and **discard** both gloves.

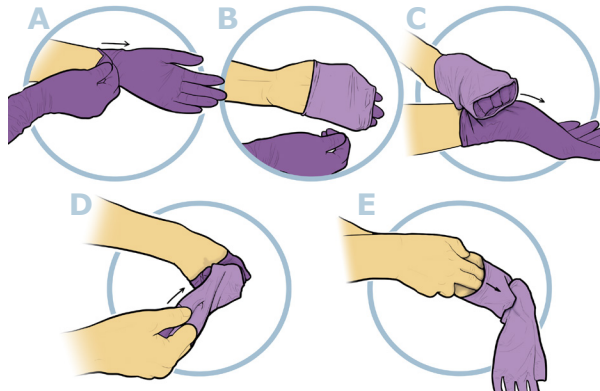


44 **Sanitize** hands and if touched, dispenser.

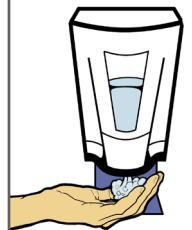


45 **Raise** arms over head in "raise the roof" motion.



46 **Observer holds** battery.

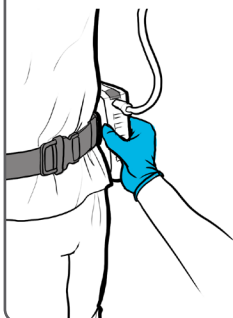

47 **Observer unbuckles** belt.

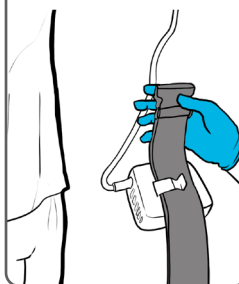

48 Observer **removes** helmet motor from behind, **avoiding** the face.



Shoe liners


49 **Put on** clean gloves




50 **Remove** shoe covers



51 **Pinch** glove cuff and **pull** over fingers, but **do not remove** glove. With covered fingers, **pinch** and **pull** opposite glove off. Using other hand, **remove** initial glove and **discard** both gloves.

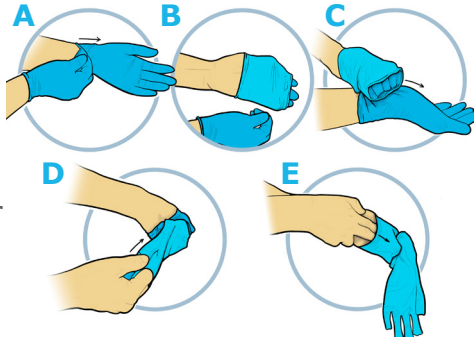

52 **Sanitize** hands and if touched, dispenser.

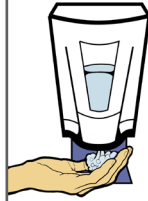

Wash

53 **Wash** hands and arms thoroughly with soap and water.



54 **Talk with your infection prevention team about showering protocols.**