

Putting PPE ON 1/1



Coveralls



PAPR system

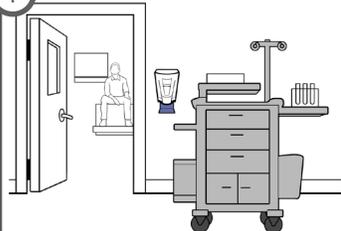


Inner gloves



Outer gloves

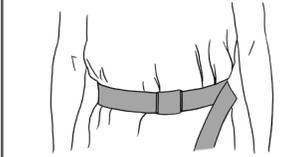
Gloves

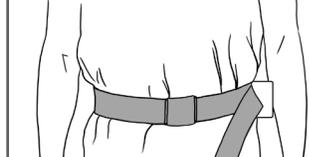
1  **Batch work and plan ahead**, having items available outside each patient room.

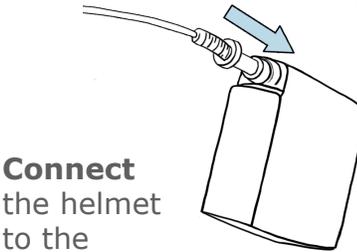
2  **Sanitize hands.**

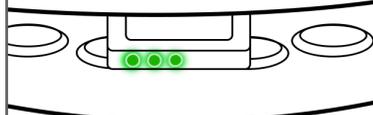
3  **Put on first pair of gloves.**

Adjust Helmet

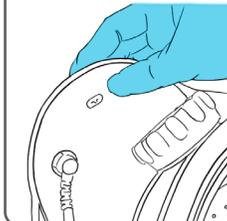
4  **Secure the belt around the waist.**

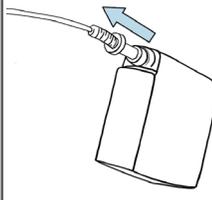
5  **Secure the battery on the left side of the belt.**

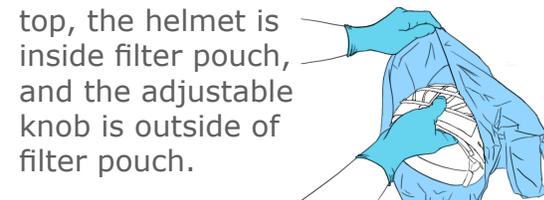
6  **Connect the helmet to the battery.**

7  **Check the LED lights to ensure the battery is fully charged.**

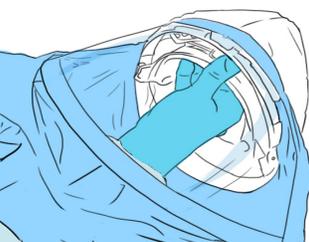
8  **Put helmet on head and size to fit by turning ratchet knob to adjust the headband.**

9  **Adjust fan speed to comfort.**

10  **Disconnect the battery.**

11  **Take off the helmet and place it inside the PAPR hood. Ensure the white HEPA filter is on top, the helmet is inside filter pouch, and the adjustable knob is outside of filter pouch.**

Put on Hood

12  **Slide the helmet into the hood until it snaps in the center hole. Snap helmet into the hood on the left and right side snaps of the face shield.**

13  **Set the assembled PAPR hood aside.**



PAPR, coveralls, double gloves

Coveralls, Gloves, PAPR

14 **Sanitize** first pair of gloves.

15 **Put on** coveralls outside room. Open-end **faces** your front.

16 **Insert** PAPR power cord connector into the battery.

17 **Put** the helmet on, letting the hood drape around your shoulders over the coveralls.

18 **Zip up** the front of the coveralls.

19 **Straighten** the hood draping and **tie** it down.

20 **Put on** second pair of gloves.

21 **Sanitize** second pair of gloves.

Entry

22 **ENTER room** **DO NOT** enter the room if you do not achieve a proper fit of the PAPR system.

23 **Do not** touch or re-adjust **PAPR** system inside the room.

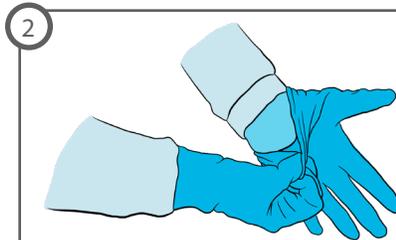


PAPR, coveralls, double gloves

Next patient: COVID+

Patient Room

1 **Sanitize** outer gloves.



Remove outer gloves.

3 **Sanitize** inner gloves.



Exit





PAPR, coveralls, double gloves

Next patient: COVID+

Gloves

After outer glove removal in the patient room, exit and...

1 **Sanitize** inner gloves.

2 **Put on** fresh outer gloves.

Walking

3 **Walk** to next patient with confirmed COVID-19 (marked) while **remaining** in marked zone.

4 **Do not touch** your face, the walls, rails, or anybody else until reaching the next room.

Entry

5 **Sanitize** outer gloves.

6 **Enter** NEXT COVID-19 patient room or area.

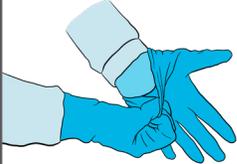
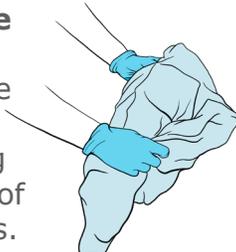
7 **Do not touch** face or re-adjust PAPR system inside room.



PAPR, coveralls, double gloves

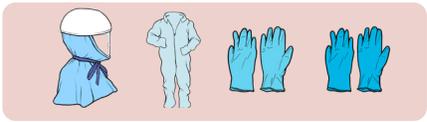
Next patient: NOT COVID+

Coveralls

- 1 Sanitize** outer gloves. 
- 2 Remove** outer gloves. 
- 3 Sanitize** inner gloves. 
- 4 Wipe down** torso, arms, and legs. 
- 5 Sanitize** inner gloves. 
- 6 Zip down** coveralls. 
- 7 Grasp** coveralls around lower back area and **pull down**. 
- 8 Remove** coveralls from shoulders and **pull down** below knees. 
- 9 Remove and step out** of the coveralls by marching and stomping on the arms of the coveralls. 
- 10 Remove** into bin for reuse without touching outside of coveralls. 
- 11 Sanitize** gloves. 

Exit

- 12 Exit** patient room. 



PAPR, coveralls, double gloves

Next patient: NOT COVID+

Hood

After gown and outer glove removal...

- Sanitize gloves.**
- Place wipe on table.**
- Untie the PAPR hood.**
- Grasp face shield at the sides and pull out snaps.**
- Pull face shield forward and release helmet pin.**
- Grab hood top, pull back, then pull forward to remove hood, avoiding your face.**
- If reusing, place hood in designated bin.**

Helmet, Belt, Battery

- Sanitize gloves.**
- Remove helmet motor, DO NOT touch your face.**
- Wipe helmet.**
- Remove and wipe battery.**
- Remove and wipe belt.**
- Sanitize gloves.**
- Remove gloves.**
- Wash hands.**

The materials are intended solely for general educational and information purposes, are made available in the context of the public health emergency related to the coronavirus (COVID-19) and have not been subject to review that typically would occur in a non-emergent situation. The materials do not constitute the provision of medical, legal or other professional advice.