

# putting PPE ON 1/2



Launderable gown



PAPR system



Inner gloves



Outer gloves

## Gloves

1 **Batch work and plan ahead**, having items available outside each patient room.

2 If you already have a mask on:

**Sanitize** hands. **Remove** mask. **Store** mask. **Sanitize** hands.

3 **Put on** first pair of gloves.

## Belt, Battery, Helmet

4 **Secure** the belt around the waist.

5 **Secure** the battery on the left side of the belt.

6 **Connect** the helmet to the battery.

7 **Check** the LED lights to ensure the battery is fully charged.

8 **Put** helmet on head and **size** to fit by turning ratchet knob to adjust the headband.

9 **Adjust** fan speed to comfort.

10 **Disconnect** the battery.

11 **Take off** the helmet and **place** it inside the PAPR hood. **Ensure** the white HEPA filter is on top, the helmet is inside filter pouch, and the adjustable knob is outside of the filter pouch.

## Hood

12 **Slide** the helmet into the hood until it snaps in the center hole. **Snap** helmet into the hood on the left and right side snaps of the face shield.

13 **Set** the assembled PAPR hood **aside**.



**PAPR, launderable gown, double gloves**

**Gown, Gloves, PAPR**

14 **Sanitize** first pair of gloves.

15 **Put on** gown outside room.  
Open-end faces your back.  
**Tie** gown at the neck.

16 **Insert** PAPR power cord connector into the battery.

17 **Tie** gown at waist.

18 **Put** the helmet on, letting the hood drape around your shoulders over the gown.

19 **Straighten** the hood draping and **tie** it down.

20 **Pull** cuffs half way over palm of the inner glove.

21 **Put on** second pair of gloves.

22 **Sanitize** second pair of gloves.

**Entry**

23 **ENTER room** **DO NOT enter** the room if you do not achieve a proper fit of the PAPR system.

24 **Do not touch or re-adjust PAPR** system inside the room.

extended  
wear **EW** <sup>1/2</sup>



**INSIDE**  
room

**PAPR, launderable gown, double gloves**

**Next patient: COVID+**

## Gloves

1 **Sanitize** outer gloves.

A hand wearing a blue nitrile glove is being held under a white hand sanitizer dispenser. The dispenser has a blue liquid inside and a nozzle at the bottom.

2

A hand wearing a blue nitrile glove is being pulled away from the wrist. The glove is being held by the other hand, which is also wearing a blue nitrile glove.

**Remove** outer gloves.

## Exit

3

**Exit** patient room.

A person wearing a PAPR and a gown is standing in a patient room. The person is walking towards a door that has an 'EXIT' sign above it. A hand sanitizer dispenser is visible on the wall next to the door.



**PAPR, launderable gown, double gloves**

**Next patient: COVID+**

**Gloves**

After outer glove removal in the patient room, exit and...

1 **Pull** cuffs half way over palm of the inner glove.

An illustration showing a person's hands. The right hand is pulling the cuff of a blue outer glove over the palm of the left hand, which is wearing a blue inner glove. An arrow points to the right, indicating the direction of the pull.

2 **Put on** fresh outer gloves.

An illustration showing a person's hands. The right hand is pulling a new blue outer glove onto the left hand, which is already wearing a blue inner glove.

**Walking**

3 **Walk** to next patient with confirmed COVID-19 (marked) while **remaining** in marked zone.

An illustration of a person in full PPE walking through a hallway. A dashed line on the floor indicates a marked zone. An arrow points forward, showing the direction of travel.

4 **Do not touch** your face, the walls, rails, or anybody else until reaching the next room.

An illustration of a person in PPE walking. Three circular icons with red 'X' marks over them show prohibited actions: touching the face, touching a wall or rail, and touching another person.

**Entry**

5 **Sanitize** outer gloves.


An illustration of a hand in a blue outer glove being placed into a white hand sanitizer dispenser.

6 **Enter** NEXT COVID-19 patient room or area.

An illustration of a person in PPE entering a room. A sign above the door says "ENTER". Inside the room, another person is sitting at a desk. A hand sanitizer dispenser is visible on the wall.

7 **Do not** touch face or re-adjust PAPR system inside room.

An illustration of a person in PPE. Red 'X' marks are placed over the face and the PAPR hood, indicating that touching the face or adjusting the hood is prohibited.

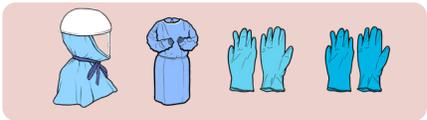

PAPR, launderable gown, double gloves

Next patient: **NOT COVID+**

## Gown, Gloves

1 **Sanitize** outer gloves.

2 **Remove** outer gloves.

3 **Sanitize** inner gloves.

4  With your back to the exit, **have a buddy stand outside the room and do the next two steps** while wearing clean gloves.

5  **Reach back and release** the neck portion.

6  **Untie** the waist.


7 **Remove** the gown into the designated bin.

8 **Sanitize** inner gloves.

## Exit

9 **Exit** patient room.

The materials are intended solely for general educational and information purposes, are made available in the context of the public health emergency related to the coronavirus (COVID-19) and have not been subject to review that typically would occur in a non-emergent situation. The materials do not constitute the provision of medical, legal or other professional advice.

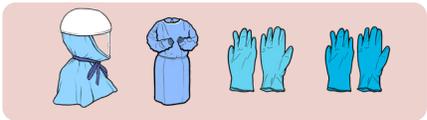

PAPR, launderable gown, double gloves

Next patient: NOT COVID+

### Hood

After gown and outer glove removal...

- Sanitize gloves.**
- Place wipe on table.**
- Untie the PAPR hood.**
- Grasp face shield at the sides and pull out snaps.**
- Pull face shield forward and release helmet pin.**
- Grab hood top, pull back, then pull forward to remove hood, avoiding your face.**
- If reusing, place hood in designated bin.**

### Helmet, Battery, Belt

- Sanitize gloves.**
- Remove helmet motor. DO NOT touch your face.**
- Put on mask.**
- Wipe helmet.**
- Remove and wipe battery.**
- Remove and wipe belt.**
- Sanitize gloves.**
- Remove gloves.**
- Wash hands.**

The materials are intended solely for general educational and information purposes, are made available in the context of the public health emergency related to the coronavirus (COVID-19) and have not been subject to review that typically would occur in a non-emergent situation. The materials do not constitute the provision of medical, legal or other professional advice.