

putting PPE ON 1/1



Coveralls



CAPR system

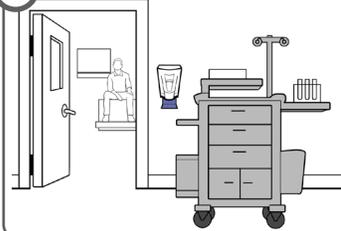


Inner gloves



Outer gloves

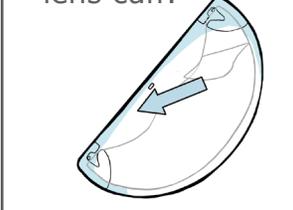
Gloves

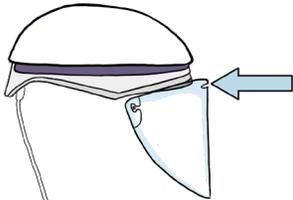
1  **Batch work and plan ahead**, having items available outside each patient room.

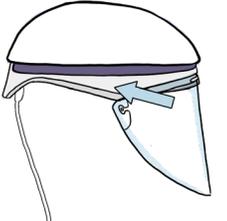
2  **Sanitize hands.**

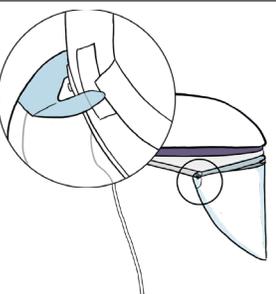
3  **Put on first pair of gloves.**

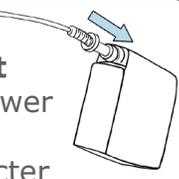
Lens Cuff

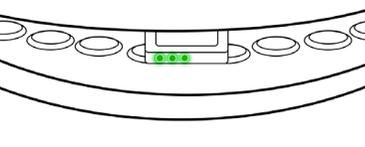
4  **Remove** the lens protector liner from the front side of the lens cuff.

5  **Secure** the lens cuff by **aligning** the center hole over the center attachment post on the helmet.

6  **Attach** the left and right side lens cuff mounting holes over the left and right side attachment posts of the helmet.
Tip: The lens cuff will click when it is properly secured.

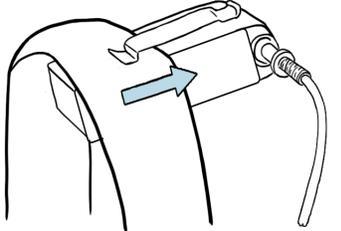
7  **Lift** the cuff flappers to rest against the foam pads at the ends of the comfort strip.

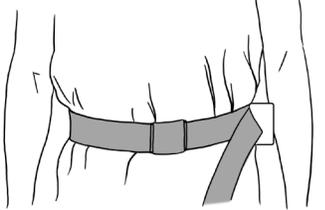
8  **Insert** the power cord connector into the battery.

9  **Check** the LED lights to ensure the battery is fully charged.

10  **Adjust** fan speed to comfort.

Battery and Belt

11  **Secure** the battery to the belt using the battery clip.

12  **Secure** the belt around the waist.



CAPR, coveralls, double gloves

Helmet

13 **Turn** ratchet knob counterclockwise to loosen the headband.

14 **Hold** the helmet near the rear headband with one hand, pull down on top of the lens cuff with the other.

15 **Place** the chin into the lens cuff while pulling the helmet over and down on the head.

16 **Turn** the ratchet knob counter clockwise to tighten the headband.

17 **Slide** the fingers around the cuff under the chin to **ensure** slight tension.

18 **Check** for proper fit, ensure:

- The cuff flappers are perpendicular to the head.
- The bottom of the headband rests 1/2 inch from eyebrows.
- LED indicators are in your peripheral vision.
- The cuff flappers are 1/4 inch from the temples.

19 **Tip:** If the headband is not in an appropriate height, doff, adjust the headband height position, and redo.

Entry

20 **Put on** coveralls outside room. Open-end **faces** your front. **Zip up** the front of the coveralls.

21 **Sanitize** first pair of gloves.

22 **Put on** second pair of gloves.

23 **Sanitize** second pair of gloves.

24 **ENTER room**
DO NOT enter the room if you do not achieve a proper fit of the CAPR system.

25 **Do not touch or re-adjust** CAPR system inside the room.

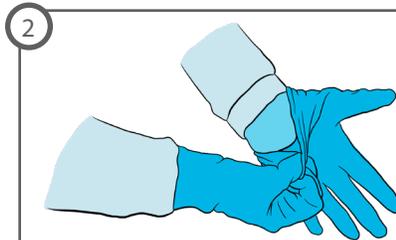


CAPR, coveralls, double gloves

Next patient: COVID+

Patient Room

1 **Sanitize** outer gloves.



Remove outer gloves.

3 **Sanitize** inner gloves.



Exit

4 **Exit** patient room.





CAPR, coveralls, double gloves

Next patient: COVID+

Gloves

After outer glove removal in the patient room, exit and...

1 **Sanitize** inner gloves.

A pair of hands wearing blue gloves is being held under a white hand sanitizer dispenser.

2 **Put on** fresh outer gloves.

A hand is shown pulling a new blue glove onto the other hand.

Walking

3 **Walk** to next patient with confirmed COVID-19 (marked) while **remaining** in marked zone.

A person in full PPE is walking in a hallway. A dashed line on the floor indicates a marked zone between two rooms.

4 **Do not touch** your face, the walls, rails, or anybody else until reaching the next room.

A person in PPE is walking in a hallway. Red 'X' marks are placed over their face, hands, and the walls to indicate that touching these areas is prohibited.

Entry

5 **Sanitize** outer gloves.

A pair of hands wearing blue gloves is being held under a white hand sanitizer dispenser.

6 **Enter** NEXT COVID-19 patient room or area.

A person in PPE is walking through a doorway into a room. A sign above the doorway says "ENTER".

7 **Do not touch** face or re-adjust CAPR system inside room.

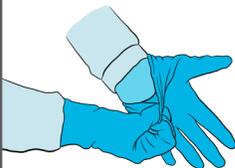
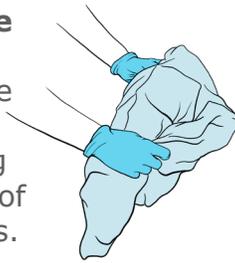
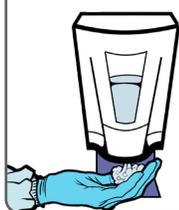
A person in PPE is shown with red 'X' marks over their face and hands, indicating that touching the face or adjusting the CAPR system is prohibited.



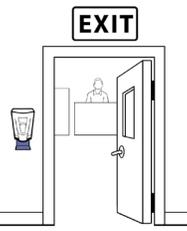
CAPR, coveralls, double gloves

Next patient: NOT COVID+

Coveralls

- 1 Sanitize outer gloves.** 
- 2 Remove outer gloves.** 
- 3 Sanitize inner gloves.** 
- 4 Wipe down torso, arms, and legs.** 
- 5 Sanitize inner gloves.** 
- 6 Zip down coveralls.** 
- 7 Grasp coveralls around lower back area and pull down.** 
- 8 Remove coveralls from shoulders and pull down below knees.** 
- 9 Remove and step out of the coveralls by marching and stomping on the arms of the coveralls.** 
- 10 Remove into bin for reuse without touching outside of coveralls.** 
- 11 Sanitize gloves.** 

Exit

- 12 Exit patient room.** 

The materials are intended solely for general educational and information purposes, are made available in the context of the public health emergency related to the coronavirus (COVID-19) and have not been subject to review that typically would occur in a non-emergent situation. The materials do not constitute the provision of medical, legal or other professional advice.



CAPR, coveralls, double gloves

Next patient: NOT COVID+

Lens Cuff

After coveralls and outer glove removal...

1 **Sanitize** gloves.

2 **Place** first wipe on table.

3 **Grab** lens cuff at the left and right side attachments. **Remove** and by pulling outward and forward.

4 **Place** lens cuff on wipe.

5 **Wipe** front and back of lens cuff. **Allow** to dry.

6 **Store** lens cuff in a gallon sized bag labeled with your name.

7 **Sanitize** gloves.

8 **Place** second wipe on table.

9 Grasp outer-shell of helmet and **take off**. **Ensure** that the helmet cord doesn't touch your body.

10 **Place** helmet on first wipe.

11 **Sanitize** gloves.

12 **Wipe** inside of the helmet.

13 **Sanitize** gloves.

14 Hold inside of the helmet and **flip** onto second wipe.

15 **Wipe** outside of the helmet.

Helmet



CAPR, coveralls, double gloves

Next patient: **NOT COVID+**

Battery

16

Sanitize gloves.

17

Wipe helmet cord.

A person in a white uniform and blue gloves is standing at a table. They are using a white cloth to wipe the cord of a helmet. A can of disinfectant is on the table.

18

Sanitize gloves.

19

Remove and wipe battery.

A person in a white uniform and blue gloves is standing at a table. They are using a white cloth to wipe a battery. A can of disinfectant is on the table.

20

Sanitize gloves.

21

Remove belt and wipe.

A person in a white uniform and blue gloves is standing at a table. They are using a white cloth to wipe a belt. A can of disinfectant is on the table.

Hands

22

Sanitize gloves.

23

Remove gloves.

A pair of hands is shown pulling apart a pair of blue nitrile gloves.

24

Wash hands.

A pair of hands is being washed under a running faucet. There is soap on the hands and bubbles in the water.